



# Captair Pyramid

Portable glove bag

Safe and convenient transport or storage of samples/evidence



## Portable glove bag

The Pyramid portable glove bag is designed to protect the user during the collection, transportation, investigation, and storage of samples or evidence.

## Easy to set up and use in the field or lab

It is a multi-function disposable glove bag made of high-quality transparent PVC. Assembled in a few seconds, it can be used anywhere. The slanted shape of the enclosure provides an ergonomic working position for the user.

Light, mobile and disposable, it is an ideal protection tool which can suit many requirements in the laboratory or in the field

## Optimal air-tightness

Each Pyramid is factory-tested to ensure an airtight seal by inflating with compressed air at a pressure of 2.5 mm Hg. Once tested, it is given a serial number and a certificate to ensure safety.



## Easy transport


Folds into an easy-to-carry size



## Ready-to-use:

- Open suspicious packages
- Collect samples on location
- Work in an inert atmosphere
- Reveal fingerprints
- Protection from splattering during biopsies
- Perform tasks while sheltered from dust and moisture
- Safely test fentanyl and unknown powders

# Specifications

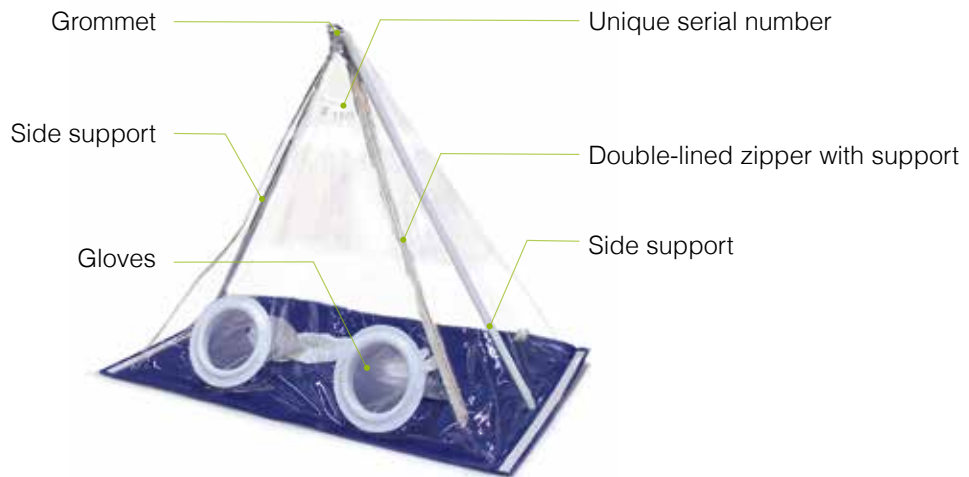

## Dimensions

Width (mm/inches)	Depth (mm/inches)	Height (mm/inches)
860 (33 <sup>7/8</sup> )	560 (22)	725 (28 <sup>1/2</sup> )

## Technical specifications

<b>Enclosure and base</b>	Flexible PVC
<b>Fastener</b>	Double lip zipper (600 mm / 23 <sup>5/8</sup> in)
<b>2 arm openings</b>	Rigid PVC (115 mm / 4 <sup>1/2</sup> in Ø)
<b>Grommet</b>	16 mm / 5/8 in Ø
<b>Operating temperature</b>	-25°C +45°C / -13°F +113°F
<b>Weight</b>	1.55 kg / 3.42 lbs

## Features

<b>Gloves</b>	Chemical Resistant Butyl rubber with PVC sleeves
<b>Serial number</b>	Unique reference number for better traceability
<b>Valve</b>	Gas insertion valve

Captair Pyramid V.03\_08\_2018 - Erlab France - Technical data are not contractual. Erlab reserves the right to modify the present document without prior notice.

France  
+33 (0) 2 32 09 55 80 | ventes@erlab.net

United States  
+1 800-964-4434 | captairsales@erlab.com

China  
+86 (0) 512 5781 4085 | sales.china@erlab.com.cn

Malaysia  
+60 (0) 7 3 555 724 | erlab@tm.net.my

Germany  
0800 330 47 31 | verkauf@erlab.net

United Kingdom  
+44 (0) 1722 341 940 | salesuk@erlab.net

Italy  
+39 (0) 2 89 00 771 | vendite@erlab.net

Spain  
+34 93 673 24 74 | ventas@erlab.net



www.erlab.com



2018 Celebrating 50 years  
of innovation in filtration