

Halo

Laboratory air filtration system*



Web user
interface access



Ceiling mounted laboratory-grade filtration unit

*Patent pending

- Removes fugitive contaminants/odors in your lab
- 24/7 filtration
- Laboratory-grade filtration using Neutrodine technology
- Low cost – no ductwork required
- Very low energy consumption (50W)
- Easy set up
- Simple operation

Healthier air for your laboratory

On the job 24/7

Halo is continuously pulling the dirty air in the lab through the filter, adsorbing the molecules that cause noxious odors, and returning clean air back to the room through low-disruptive side vents. The laboratory-grade molecular filtration system contains Neutrodine, a unique and powerful activated carbon mixture combined with a pre-filter that captures a large variety of airborne contaminants, which left untreated, could represent a potential health hazard.


Simple to install

Halo is easily installed in the ceiling - no need to connect to your HVAC. Multiple units can be installed to increase the coverage area.

Easy to operate


A soft band of light assures the occupants of the lab that the air is clean. When the LED lights pulse, Halo is letting you know that it is sensing and capturing fugitive contaminants. You can easily check the status, change the setting, or check for routine maintenance via your smartphone.



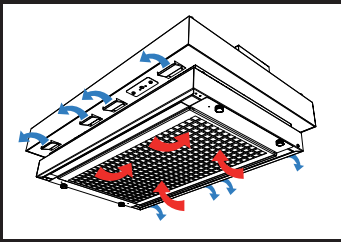

• Fugitive contaminants are pulled into the filter, molecules are captured, and clean air is recirculated back into the room.

Dimensions

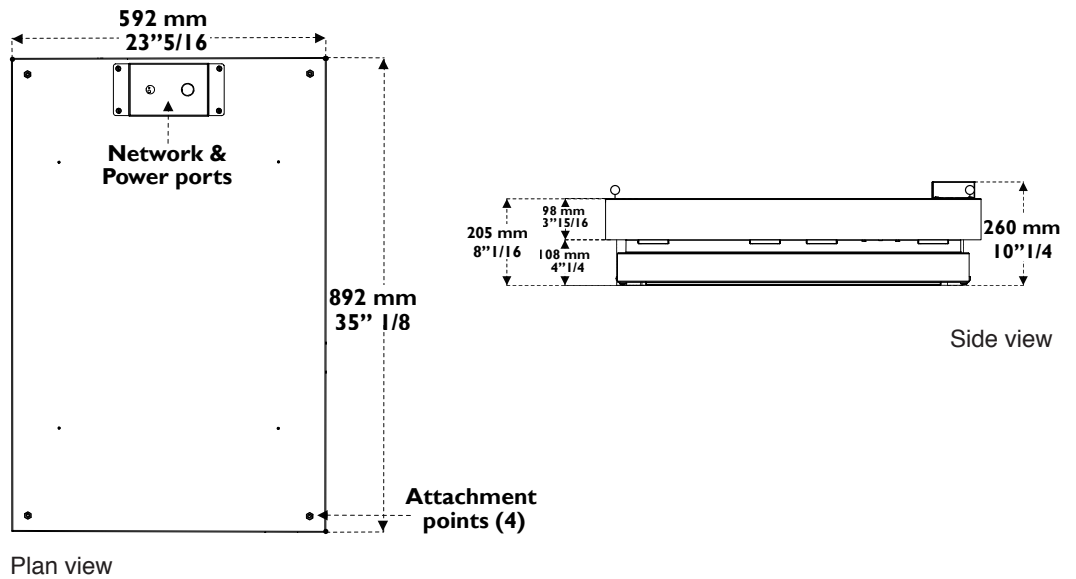

Specifications

CFM	130 ft ³ - 220 m ³ /h
Coverage area	250 square feet
Operating modes	24/7, Night/Day, Alert
Pre filtration	Particulate prefilter
Filter options (using Neutrodine Filtration technology)	- VOC : for volatile organic compounds - Chemplus : for broad range of chemicals
e-Guard Communication	- Webservice (Embedded) - Software (Start & Extended versions) - BACnet
Ports	1 x Ethernet Port
Installation	Hung via 4 eye bolts (included)
Weight	68 lbs - 31 kg (including filter)
Voltage / frequency	80-240V 50/60hz
Energy consumption	50 W
Decibels	Minimum flow : 45 / Maximum flow : 60

About Erlab

Inventor of the ductless laboratory fume hood and world-wide leader since 1968, Erlab's passion is to focus on the research & development, design, and manufacturing of cutting-edge toxic gas air filtration in the laboratory.

As an innovator, Erlab is committed to safety, performance, energy efficiency and sustainability and has remained number one in the world for ductless fume hoods since 1968 with more than 100 000 units in operation.

www.erlab.com

North America
+1 800-964-4434 | captairsales@erlab.com

China
+86 (0) 512 5781 4085 | sales.china@erlab.com.cn

France
+33 (0) 2 32 09 55 80 | ventes@erlab.net

Germany
0800 330 47 31 | verkauf@erlab.net

Italy
+39 (0) 2 89 00 771 | vendite@erlab.net

Malaysia
+60 (0) 7 3 555 724 | erlab@tm.net.my

Spain
+34 93 673 24 74 | ventas@erlab.net

United Kingdom
+44 (0) 1722 341 940 | salesuk@erlab.net