

## Product datasheet

# Captair 834 Smart & Midcap

Ductless filtering chemical storage cabinets

### Safer to operate

- Exclusive Erlab filtration technology
- Filtration performances tested in accordance with AFNOR NF X 15 211
- Secured storage
- Real time sensors to detect main filter saturation with solvents, acids or formaldehyde\*
- Erlab Safety Program: be guided for the right choice of your secured storage solutions
- Connect with eGuard™ for real-time status, notifications and alerts\*

### Simpler to use with Smart-Light™

With Smart-Light, you can easily see that the hood is operating safely.

Should the light pulse, you are notified that\*:

- Containment has been compromised or,
- The filter has breakthrough or,
- There is a fan failure

### Flexibility

- 2 possible configurations:  
Double doors and shelves  
Pull-out doors with storage trays

### Savings

- No ductwork needed
- Low annual energy cost
- Low replacement filter cost

### Environment

- Pollutants filtered at the source

\* Not available on Captair 834 Midcap



You get the highest level of filtration performance



Smart-Light keeps you safe





Connect your equipment and download the eGuard App for remote access




Download our eGuard application





Option 1 – Swing doors with shelves

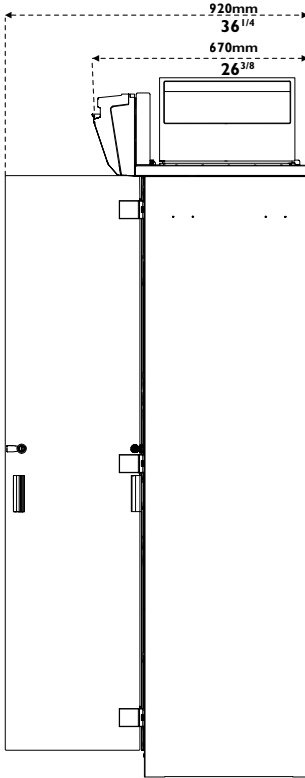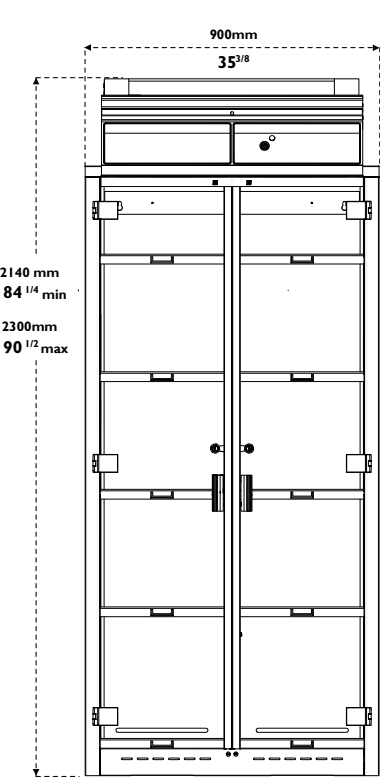


Shelf with built-in spill retention tray

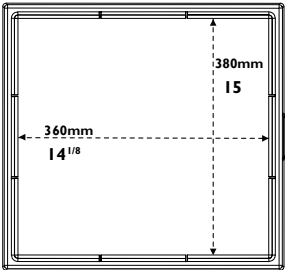

Option 2 – Pull-out doors with storage trays

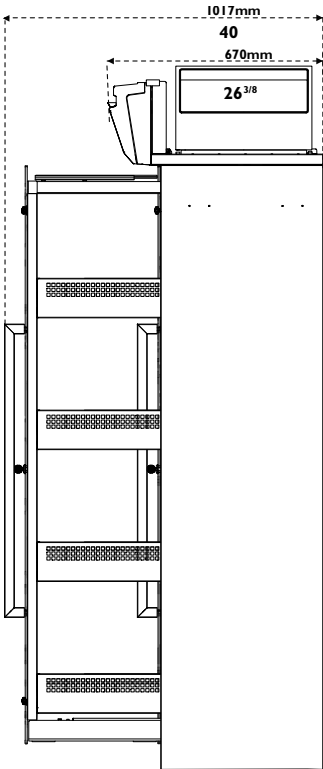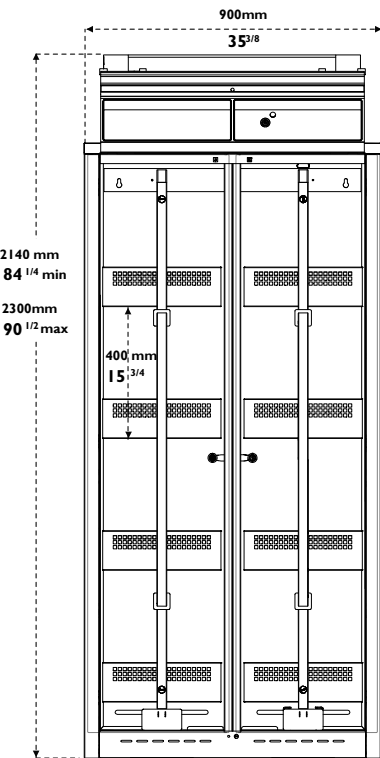





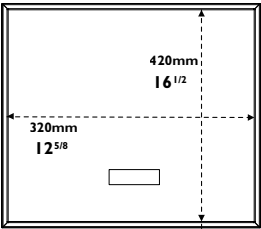



Storage with built-in spill retention tray





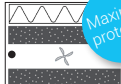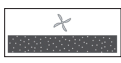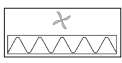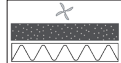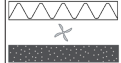

Designed with you in mind: our filtration column can be configured for your specific application requirements

		Products handled / Applications			
		Liquid chemicals handlings	Powders handlings	Liquid chemicals and powders handlings	Liquid chemicals handlings in clean room
Customized filtration column	Class 1 according to the NF X 15-211	 2C	NA	 1P2C	 2C1P
	Class 2 according to the NF X 15-211	 1C	 1P	 1P1C	 1C1P

**C**  **Carbon filtration for gases and vapours**  
AS:For organic vapours  
BE+:Polyvalent for acid + organic vapours  
F:For formaldehyde vapours  
K:For ammonia vapours

**P**  **Particulate filtration for powders**  
HEPA H14:99.995 % efficiency filtration of particles over 0.1 µm in size  
ULPA U17:99.999995 % efficiency filtration of particles over 0.1 µm in size

 **Ventilation**

• **Molecode**  
Automatic alarm to detect filter saturation

Safety Standards	Filtration performances tested according to the AFNOR NF X 15-211:2009 standard : France EN 1822 : 1998 (HEPA H14 & ULPA U17 Filters) - CE Marking
Air Flow	220 m3/h / 129 CFM
Voltage/Frequency	80 - 220 V / 50-60 Hz
Power consumption	45 W
Structure	Corrosion resistant electro-galvanized steel coated with anti-acid polymer
Doors	Clear, chemical resistant acrylic for easy viewing
Filtration Module	Polypropylene

## Features

Smart-Light*	When the light is pulsing: -Door(s) left open -Containment is compromised -Filter breakthrough
Filtration technology	1 column that can be configured to handle liquids, powders, or both
Carbon filtration for gases and vapors	Depending on the filtration column configuration (see above)
Particulate filtration for powders	Depending on the filtration column configuration (see above)
eGuard app* (Android or iOS)	For remote control to monitor the status of the storage cabinet, change the settings, and deliver safety alerts immediately to your devices (mobile, tablet and PC)
Doors sensors*	Alarm if doors are left open
Chemical Listing	List of most commonly used chemicals (NIOSH) and their retention capacity by our filters

## Storage configurations

	Option 1 – Double doors with shelves	Option 2 – Pull-out doors with storage trays
Storage capacities	120 x 1L bottles	100 x 1L bottles
Storage compartments	2	2
Delivered with	10 adjustable shelves with integrated retention tray	8 fixed trays
Absorbing mats	2	8
Lock	Key lock	

## Options

Molécode*	Detection sensor for : Type S, for solvents / Type A, for acids / Type F, for formaldehydes
-----------	---

\* Not available on Captair 834 Midcap



*Erlab's state of the art Research & Development Laboratory relies exclusively on filtration*

## About Erlab

### We provide safety, we protect your health

Erlab invented the ductless fume hood in 1968. With more than 50 years of experience in the field of chemical filtration and protection of laboratory personnel; we know the formula for safety. With Erlab, you will never have to wonder or worry if our products are safe. We build each one of the following 7 ingredients into our products, and without all of them, your health and safety will be compromised.

#### 1 Erlab R&D Laboratory

The engineers and chemists in our state-of-the-art R&D laboratory understand molecular filtration. We are committed to designing products that are safe and of the highest quality, strive to improve our products, and continuously develop new products that provide greater protection in the laboratory.

#### 2 Strict Safety Standards

We hold ourselves to the highest standard and adhere to the strict AFNOR NF X 15-211: 2009 filtration safety standard as endorsed by ANSI Z9.5-2012.

#### 3 A Published Chemical Listing

It all begins here. Without this listing, we are not compliant with AFNOR NFX 15-211. Our in-house laboratory tests, as well as independent testing, to verify the retention capacity of over 700 chemicals for our filters.

#### 4 Independent Testing

Erlab filters have been independently tested multiple times at various concentrations guaranteeing that our safety solutions all adhere to the strict performance criteria of the AFNOR NF X 15-211:2009 standard assuring that the emission concentration at the filter exhaust will always be lower than 1% of the TLV.

#### 5 Application Questionnaire (Valiquest)

Our laboratory specialists will recommend the appropriate filtration fume hood, type of filter, and personalized advice.

#### 6 Certificate of Validation for the chemicals used in the hood

A certified PhD chemist issues a Certificate of Validation with a list of the chemicals approved for use in the hood.

#### 7 Our Safety Program

We back up our products 100%. This program includes your specialized chemical evaluation, validation of your hood upon installation, and a filtration safety specialist at your service to ensure that your hood is operating to its full potential.

France  
+33 (0) 2 32 09 55 80 | ventes@erlab.net

United States  
+1 800-964-4434 | captairsales@erlab.com

China  
+86 (0) 512 5781 4085 | sales.china@erlab.com.cn

Spain  
+34 936 732 474 | export.south@erlab.net

Germany  
0800 330 47 31 | export.north@erlab.net

United Kingdom  
+44 (0) 1722 341 940 | export.north@erlab.net

Italy  
+39 (0) 2 89 00 771 | export.south@erlab.net